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
**Counterfeit medical products seized by medical Devices
Sector inspectors**

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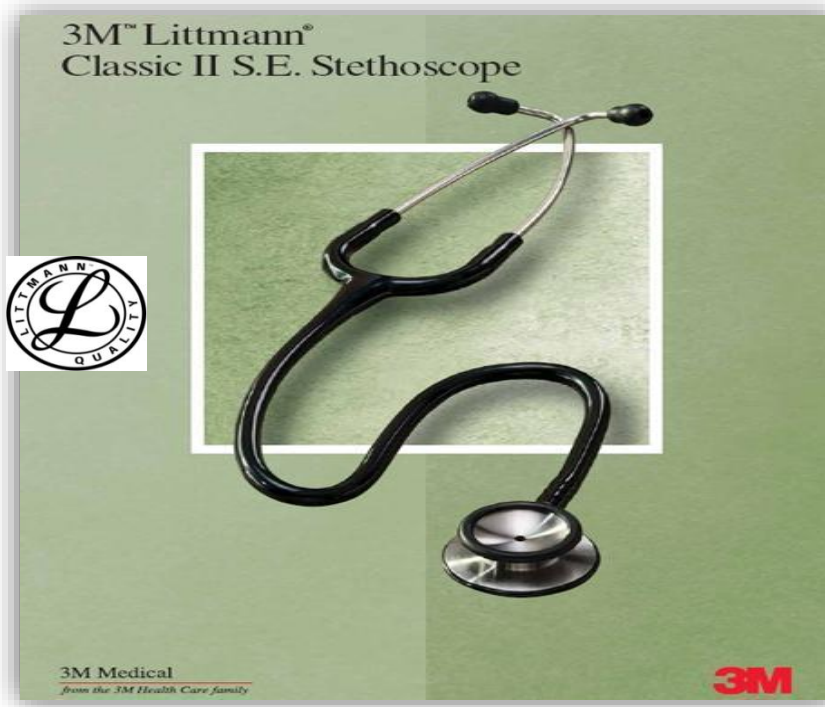Saudi Food and Drug Authority , would like to warn all hospitals , clinics , Health care practitioners of the existence of a counterfeit **Littmann classic II S.E stethoscope** marketed as the original product manufactured by **3M** Health care according to following information :-

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The SFDA has become aware through 3M scientific office for medical care in Riyadh of suspected counterfeit **Littmann classic S.E.II stethoscope** devices. Preliminary testing has shown that the counterfeit stethoscope deliver correspondingly poor performance comparing to the Stethoscope manufactured by 3M Health care and that the materials used in manufacturing it, is of inferior quality to those used in the genuine stethoscope .

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Instructions:

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The SFDA would like from all hospitals, clinics, Health care practitioners taking the following instructions:

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- Make sure that you buy directly from the manufacturer or an authorized dealer / distributor.
- Replace counterfeit with the genuine device.
- Complete the warranty form (if provided) so that the manufacturer can contact buyer if they need to (e.g. in the event of a field safety corrective action / recall).
- Report to SFDA through NCMDR website <http://ncmdr.sfda.gov.sa> if there is any suspected counterfeit device/product or call 19999
- Do not buy medical devices from internet

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Please refer to the figures below to view the differences between the genuine and the counterfeit stethoscope.

Counterfeit Medical Device جهاز طبي مقلد

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Saudi Food and Drug Authority warns all hospitals , clinics , Health care practitioners of the existence of a counterfeit Laryngoscope set (kit) devices marketed as the original product manufactured by

Riester according to following information :-

Counterfeit Medical Device جهاز طبي مقلد

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Port of entry inspectors in SFDA have detected a counterfeit Laryngoscope . Preliminary testing has shown that There are differences between the original product manufactured by **Riester** and the counterfeit device in labeling, packaging and designing of sets.

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Instructions:

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The SFDA would like from all hospitals , clinics , Health care practitioners taking the following instructions.

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- Ensure that buying directly from the manufacturer or an authorized dealer / distributor.
- Replace counterfeit with the genuine device.
- Complete the warranty form (if provided) so that the manufacturer can contact buyer if they need to (e.g. in the event of a field safety corrective action / recall).
- Report to SFDA through NCMDR website <http://ncmdr.sfda.gov.sa> if there is any suspected counterfeit device. or call 19999
- Don't buy medical devices from internet

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Please refer to the figures below to view the differences between the genuine and the counterfeit Laryngoscope Set.

Counterfeit Medical Device جهاز طبي مقلد

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Saudi Food and Drug Authority , would like to warn all hospitals , clinics , Health care practitioners of the existence of a counterfeit **Diagnostic Set** used for to look into the eyes, ears, nose, and throat, marketed as the original product manufactured by  according to following information :-

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A counterfeit diagnostic set have been detected by SFDA Port of entry inspectors which are used to look into the eyes, ears, nose, and throat . Preliminary testing has shown that there are differences between original Diagnostic set manufactured by **Riester** and counterfeit in labeling, packaging and designing of sets .

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- Ensure to buy directly from the manufacturer or an authorized dealer / distributor.
- Replace counterfeit with the genuine device.
- Complete the warranty form (if provided) so that the manufacturer can contact buyer if they need to (e.g. in the event of a field safety corrective action / recall).
- Report to SFDA through NCMDR website <http://ncmdr.sfda.gov.sa> if there is any suspected counterfeit device .
- Do not buy medical devices from internet

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Please refer to the figures below to view the differences between the genuine and the counterfeit Diagnostic Set

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Saudi Food and Drug Authority , would like to warn all hospitals, Health care practitioners of the existence of a counterfeit Endodontic Files devices marketed as the original product manufactured by **DENSPLY** according to following information :-

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The SFDA has become aware of the counterfeit **endodontic** files. The suspected counterfeit devices have subtle differences in labeling, packaging and designing of files regarding the original device manufacturer by **Dentsply** . A counterfeit endodontic files may break inside the tooth during the treatment due to poor quality.

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Instructions:

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Saudi Food and Drug Authority ,would like to warn all hospitals, Health care practitioners to follow instructions:

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- Make sure that you buy directly from the manufacturer or an authorized dealer / distributor.
- Replace counterfeit with the genuine device.
- Complete the warranty form (if provided) so that the manufacturer can contact you if they need to (e.g. in the event of a field safety corrective action / recall).
- Report to SFDA through NCMDR website <http://ncmdr.sfda.gov.sa> if there is any suspected counterfeit device or call 19999
- Do not buy medical devices from internet

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Please refer to the figures below to view the differences between the genuine and the counterfeit

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Saudi Food and Drug Authority , would like to warn all hospitals , clinics , Health care practitioners of the existence of a counterfeit **stethoscopes** marketed as the original product "**duplex deluxe stethoscope**" manufactured by **Riester** according to following information :-

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Port of entry inspectors in SFDA have detected a counterfeit **duplex deluxe** stethoscope . Preliminary testing has shown that the counterfeit duplex deluxe stethoscope deliver correspondingly poor performance comparing to the Stethoscope manufactured by **Riester** and that the materials used in manufacturing it, is of inferior quality to those used in the genuine stethoscope .

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- Ensure that buying directly from the manufacturer or an authorized dealer / distributor.
- Replace counterfeit with the genuine device.
- Complete the warranty form (if provided) so that the manufacturer can contact buyer if they need to (e.g. in the event of a field safety corrective action / recall).
- Report to SFDA through NCMDR website <http://ncmdr.sfda.gov.sa> if there is any suspected counterfeit device or call 19999
- Don't buy medical devices from internet .

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Please refer to the figures below to view the differences between the genuine and the counterfeit **Stethoscope**.

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